

# HeartSine® samaritan® PAD 450P AED

Automated external defibrillator with integrated  
CPR Rate Advisor™

## Data sheet

### Key link in the chain of survival

Cardiopulmonary Resuscitation (CPR) and Automated External Defibrillators (AEDs) are key links in the chain of survival of sudden cardiac arrest (SCA). Some cardiac events are treatable with effective CPR alone. Others require a combination of effective CPR and the delivery a lifesaving shock by an AED. Either way, every minute counts.

Only about five percent of SCA victims survive. However, survival rates can increase up to 74%<sup>1</sup> if CPR and a shock from an AED are provided within three minutes of collapse. Reducing response time by even one or two minutes from collapse to shock can mean the difference between death and survival.<sup>2</sup>

Offering real-time CPR Rate feedback, the HeartSine samaritan PAD 450P (SAM 450P) Automated External Defibrillator (AED) with integrated CPR Rate Advisor meets the needs of two key links in the chain of survival. Not only can the semi-automatic SAM 450P deliver a lifesaving shock, it provides real-time visual and verbal feedback to the rescuer on the rate of CPR compressions during an SCA resuscitation — effectively assisting the rescuer to perform CPR.



## Real-time CPR rate feedback



### Integrated real-time CPR rate feedback

Easy-to-understand visual and voice prompts guide the rescuer through the entire resuscitation process, providing specific feedback on the rate of compressions (“Push faster,” “Push slower,” “Good speed”).



### Reminders to perform CPR

To encourage more hands-on time for CPR delivery, the HeartSine samaritan PAD 450P continues to remind the rescuer to perform CPR when no CPR is detected.

## Ready to shock



### Unique Pediatric-Pak

Ensures the guidelines-recommended energy level is delivered for children, between 1 and 8 years of age or up to 55 lb (25 kg).



### High level of protection from dust and water

Offers IP56 rating, one of the highest ratings in the industry.



### Clinically validated technology<sup>3</sup>

Advanced electrode technology and SCOPE (Self-Compensating Output Pulse Envelope) biphasic technology, a low energy escalating waveform that automatically adjusts for differences in patient impedance.



### Highly portable

With the lightest weight and most compact footprint among leading AEDs, is easily transported and fit into constrained spaces.

## Simple to own



### Two parts, one expiration date

The innovative Pad-Pak, an integrated battery and electrode single-use cartridge with one expiration date, offers one simple maintenance change every four years.




### Low cost of ownership

Shelf life of four years means that the Pad-Pak may offer savings over other defibrillators that require separate battery and electrode replacements.



### 8-year warranty

AED is backed by an 8-year limited warranty.

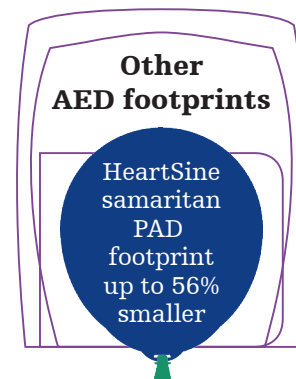

# Specifications

## Defibrillator

---

**Waveform:** Self-Compensating Output Pulse Envelope (SCOPE) optimized biphasic escalating waveform compensates energy, slope and duration for patient impedance

## Patient analysis system

---

**Method:** Evaluates patient's ECG, electrode contact integrity and patient impedance to determine if defibrillation is required

**Sensitivity/Specificity:** Meets IEC/EN 60601-2-4

**Impedance range:** 20-230 ohms

## Energy selection

---

### Pad-Pak:

Shock 1: 150 J  
Shock 2: 150 J  
Shock 3: 200 J

### Pediatric-Pak:

Shock 1: 50 J  
Shock 2: 50 J  
Shock 3: 50 J

### Charge time (typical):

150J in < 8 seconds  
200J in < 12 seconds

## Environmental

---

### Operating/Standby temperature:

32°F to 122°F (0°C to 50°C)

### Transport temperature:

32°F to 122°F (0°C to 50°C)

NOTE: It is recommended that the device should be placed in an ambient temperature of between 32°F to 122°F (0°C to 50°C) for at least 24 hours upon first receipt.

**Relative humidity:** 5% to 95% non-condensing

**Water resistance:** IEC 60529/ EN60529 IPX6 with electrodes connected and battery installed

**Dust resistance:** IEC 60529/ EN60529 IP5X with electrodes connected and battery installed

**Enclosure:** IEC/EN 60529 IP56

**Altitude:** -1,250 to 15,000 feet (-381 to 4,575 meters)

**Shock:** MIL STD 810F Method 516.5, Procedure 1 (40 G's)

**Vibration:** MIL STD 810F Method 514.5, Procedure 1

Category 4 Truck Transportation – US Highways

Category 7 Aircraft – Jet 737 & General Aviation

**Atmospheric pressure:** 572 hPa to 1060hPa (429 mmHg to 795 mmHg)

**EMC:** IEC/EN 60601-1-2

**Radiated emissions:** IEC/EN 55011

### Electrostatic discharge:

IEC/EN 61000-4-2 (8 kV)

### RF immunity:

IEC/EN 61000-4-3 80MHz-2.5 GHz, (10 V/m)

### Magnetic field immunity:

IEC/EN 61000-4-8 (3 A/m)

**Aircraft:** RTCA/DO-160G, Section 21 (Category M)

RTCA/DO-227 (TSO/ETSO-C142a)

**Falling height:** 3.3 feet (1 meter)

## Physical characteristics

---

With Pad-Pak inserted:

### Size:

8.0 in x 7.25 in x 1.9 in  
(20 cm x 18.4 cm x 4.8 cm)

**Weight:** 2.4 lb (1.1 kg)

## Accessories

---

### Pad-Pak Electrode and Battery Cartridge

**Shelf life/Standby life:** See the expiration date on the Pad-Pak/Pediatric-Pak (4 years from manufacture date)

**Weight:** 0.44 lb (0.2 kg)

### Size:

3.93 in x 5.24 in x 0.94 in  
(10 cm x 13.3 cm x 2.4 cm)

**Battery type:** Disposable single-use combined battery and defibrillation electrode cartridge (lithium manganese dioxide (LiMnO<sub>2</sub>) 18V)

### Battery capacity (new):

> 60 shocks at 200J or 6 hours of battery use

**Electrodes:** Disposable defibrillation pads are supplied as standard with each device

**Electrode placement:** Anterior - lateral (Adult)

Anterior - posterior or Anterior - lateral (Pediatric)

**Electrode active area:** 15 in<sup>2</sup> (100 cm<sup>2</sup>)

**Electrode cable length:** 3.3 feet (1 meter)

**Aircraft safety test (TSO/ETSO-certified Pad-Pak):** RTCA/DO-227 (ETSO-C142a)

## Data storage

---

**Memory type:** Internal memory

**Memory storage:** 90 minutes of ECG (full disclosure) and event/incident recording

**Review:** Custom USB data cable (optional) directly connected to PC with Saver EVO Windows-based data review software

## Materials used

---

**Defibrillator housing:** ABS, Santoprene

**Electrodes:** Hydrogel, Silver, Aluminium and Polyester

## Warranty

---

**AED:** 8-year limited warranty

## References

1. Valenzuela TD, et al. 2000. Outcomes of Rapid Defibrillation by Security Officers After Cardiac Arrest in Casinos. *New England Journal of Medicine*. 343:1206-09.
2. Mosesso Jr VN, et al. 2002. Proceedings of the National Center for Early Defibrillation Police AED Issues Forum. *Prehospital Emergency Care*. 6(3):273-82.
3. Walsh SJ, McClelland A, Owens CG, et al. Efficacy of distinct energy delivery protocols comparing two biphasic defibrillators for cardiac arrest. *Am J Cardiol*. 2004;94:378-380.

## HeartSine® samaritan® PAD Automated External Defibrillators (AEDs)

### BRIEF SUMMARY OF INDICATIONS AND IMPORTANT SAFETY INFORMATION


**INDICATIONS FOR USE:** The HeartSine samaritan PAD SAM 350P (SAM 350P), HeartSine samaritan PAD SAM 360P (SAM 360P) and HeartSine samaritan PAD SAM 450P (SAM 450P) are indicated for use on victims of cardiac arrest who are exhibiting the following signs: unconscious, not breathing, without circulation (without a pulse). The devices are intended for use by personnel who have been trained in their operation. Users should have received training in basic life support/AED, advanced life support or a physician-authorized emergency medical response training program. The devices are indicated for use on patients greater than 8 years old or over 55 lb (25 kg) when used with the adult Pad-Pak (Pad-Pak-01 or Pad-Pak-07). They are indicated for use on children between 1 and 8 years of age or up to 55 lb (25 kg) when used with the Pediatric-Pak (Pad-Pak-02).

**CONTRAINDICATION:** If the patient is responsive or conscious, do not use the HeartSine samaritan PAD to provide treatment.

**WARNINGS: AEDs:**

- The HeartSine samaritan PAD delivers therapeutic electrical shocks that can cause serious harm to either users or bystanders. Take care to ensure that no one touches the patient when a shock is to be delivered.
- Touching the patient during the analysis phase of treatment can cause interference with the diagnostic process. Avoid contact with the patient while the HeartSine samaritan PAD is analyzing the patient. The device will instruct you when it is safe to touch the patient.
- Do not delay treatment trying to find out the patient's exact age and weight. If a Pediatric-Pak or an alternative suitable defibrillator is not available, you may use an adult Pad-Pak.
- The SAM 360P is a fully automatic defibrillator. When required, it will deliver a shock to the patient WITHOUT user intervention.
- The SAM 450P CPR Rate Advisor is currently only intended to provide feedback on adult patients. If you treat a pediatric patient with the SAM 450P and an adult Pad-Pak, ignore any voice prompts regarding the rate of CPR.
- Do NOT use the HeartSine samaritan PAD in the vicinity of explosive gases, including flammable anesthetics or concentrated oxygen.
- Do NOT open or repair the device under any circumstances as there could be danger of electric shock. If damage is suspected, immediately replace the HeartSine samaritan PAD.

**Pad-Paks:**

- Do not use if the gel is dry.
- The Pediatric Pad-Pak is not for use on patients under 1 year old. For use with children up to the age of 8 years or up to 55 lb (25 kg). DO NOT DELAY THERAPY IF YOU ARE NOT SURE OF EXACT AGE OR WEIGHT.
- Only HeartSine samaritan PADs with the  label are suitable for use with the Pediatric-Pak. If the HeartSine samaritan PAD you are using does not have this label, use the adult Pad-Pak if no alternatives are available.
- The use of the Pediatric-Pak will enable delivery of 50J shocks to the pediatric patient.
- The Pediatric-Pak contains a magnetic component (surface strength 6500 gauss). Avoid storage next to magnetically sensitive storage media. It is advised that Pediatric-Paks are stored separately when not in use.
- Never charge, short circuit, puncture, deform, incinerate, heat above 85°C or expose contents of TSO (Aviation) Pad-Pak to water. Remove when discharged.

**PRECAUTIONS: AEDs:**

- Proper placement of the HeartSine samaritan PAD electrode pads is critical. Electrode pads must be at least 1 in (2.5 cm) apart and should never touch one another.
- Do not use electrode pads if pouch is not sealed.
- Check the device periodically in accordance with the service and maintenance instructions provided in the User Manual.
- Operate the HeartSine samaritan PAD at least 6 feet (2 meters) away from all radio frequency devices or switch off any equipment causing interference.
- Use of the device outside the operating and storage ranges specified in the User Manual may cause the device to malfunction or reduce the shelf life of the Pad-Pak.
- Do not immerse any part of the HeartSine samaritan PAD in water or any type of fluid.
- Do not turn on the device unnecessarily as this may reduce the standby life of the device.
- Do not use any unauthorized accessories with the device as the HeartSine samaritan PAD may malfunction if non-approved accessories are used.
- Dispose of the device in accordance with national or local regulations.
- Check with the relevant local government health department for information about any requirements associated with ownership and use of a defibrillator in the region where it is to be used.

**Pad-Paks:**

- Check expiration date.

**Saver EVO Software:**

- Download the complete HeartSine samaritan PAD memory prior to erasing it. This information should be stored safely for future reference. Ensure that only the events you want to delete have been selected prior to deleting. Once deleted from your computer's memory, events cannot be regenerated and all information will be lost.

**POTENTIAL ADVERSE EFFECTS:** The potential adverse effects (e.g., complications) associated with the use of an automated external defibrillator include, but are not limited to, the following:

- Failure to identify shockable arrhythmia.
- Failure to deliver a defibrillation shock in the presence of VF or pulseless VT, which may result in death or permanent injury.
- Inappropriate energy which could cause failed defibrillation or post-shock dysfunction.
- Myocardial damage.
- Fire hazard in the presence of high oxygen concentration or flammable anesthetic agents.
- Incorrectly shocking a pulse-sustaining rhythm and inducing VF or cardiac arrest.
- Bystander shock from patient contact during defibrillation shock.
- Interaction with pacemakers.
- Skin burns around the electrode placement area.
- Allergic dermatitis due to sensitivity to materials used in electrode construction.
- Minor skin rash.

**CAUTION:** U.S. Federal law restricts this device to sale by or on the order of a physician.

**Please consult the user manual at [heartsine.com](http://heartsine.com) for the complete list of indications, contraindications, warnings, precautions, potential adverse events, safety and effectiveness data, instructions for use and other important information.**

All claims valid as July 2021.

**For further information, please contact Stryker at 800 442 1142 (U.S.) or visit our website at [strykeremergencycare.com](http://strykeremergencycare.com)**

## Emergency Care Public Access

Stryker's AEDs require a prescription in the U.S. Please consult your physician. AED users should be trained in CPR and in the use of the AED.

Although not everyone can be saved, studies show that early defibrillation can dramatically improve survival rates. AEDs are indicated for use on adults and children. AEDs may be used on children weighing less than 55 lb (25 kg) but some models require separate defibrillation electrodes.

The information presented is intended to demonstrate Stryker's product offerings. Refer to operating instructions for complete directions for use indications, contraindications, warnings, cautions, and potential adverse events, before using any of Stryker's products. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your representative if you have questions about the availability of Stryker's products in your area. Specifications subject to change without notice.

Stryker or its affiliated entities own, use, or have applied for the following trademarks or service marks: CPR Rate Advisor, HeartSine, Pad-Pak, Pediatric-Pak, samaritan, Saver EVO, SCOPE, Stryker. All other trademarks are trademarks of their respective owners or holders.

The absence of a product, feature, or service name, or logo from this list does not constitute a waiver of Stryker's trademark or other intellectual property rights concerning that name or logo.

 Pad-Pak; Pediatric-Pak

 HeartSine samaritan PAD: UL Classified. See complete marking on product.

Date of Issue: 07/2021; Made in U.K.; H009-041-001-AG EN-US  
HeartSine SAM 450P is not available for sale outside of the U.S. or Japan.  
Copyright © 2021 Stryker.



### Manufactured by:

HeartSine Technologies Ltd.  
207 Airport Road West  
Belfast, BT3 9ED  
Northern Ireland  
United Kingdom  
Tel +44 28 9093 9400  
Fax +44 28 9093 9401  
[heartsinesupport@stryker.com](mailto:heartsinesupport@stryker.com)  
[heartsine.com](http://heartsine.com)

### Distributed in U.S. by:

Stryker Emergency Care  
11811 Willows Road NE  
Redmond, WA, 98052 U.S.A.  
Toll free 800 442 1142  
[strykeremergencycare.com](http://strykeremergencycare.com)