



RADIAN QBio Co., Ltd.

www.radianqbio.ae

HEART GUARDIAN HR-701 PLUS


WITH MONITOR

PATIENT REAL TIME ECG & VIDEO PROMPTS



AED Defibrillators Machine **Dubai**
Dubai - UAE
E-mail: info@radianqbio.ae
Website: <http://www.radianqbio.ae>

HR-701 Plus

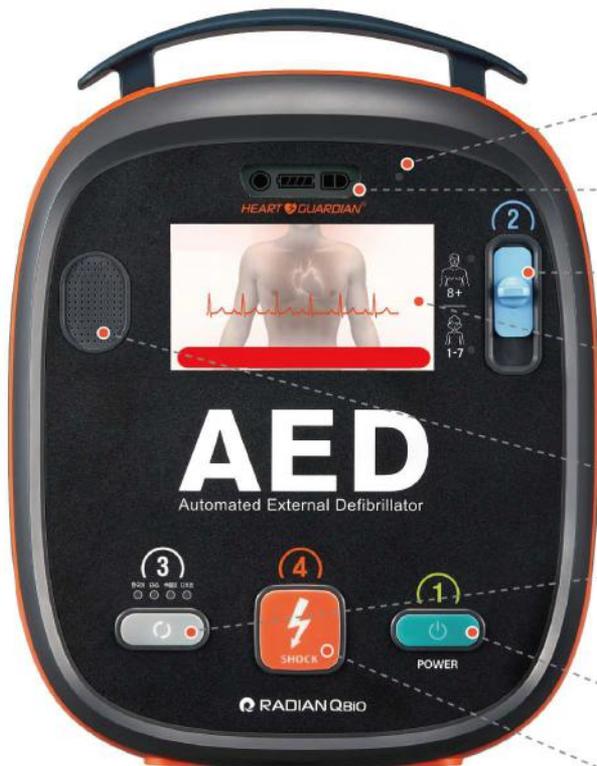


▲ HR-701 Plus

- ▶ Ambient Light Sensor
- ▶ Status Indicators (Operation / Battery / Pads)
- ▶ Adult / Child Mode by Slide Switch
- ▶ 4.3" LCD display
(Video Prompts / Patient Real Time ECG)
- ▶ Speaker (Voice Prompts)
- ▶ Language Select button (Customized Languages)
- ▶ Power Button
- ▶ Shock Button




AED HR-701 Plus

Heart Guardian

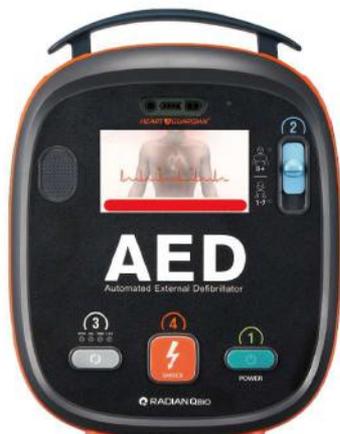

- Automatic Self-test once a month
- Adult / Child Mode Compatible switch (Energy: 150J~250J for Adult / 50J for Child)
- 4.3" High Resolution Color LCD display
- Size: 222(W) x 280(D) x 75(H)mm
- Weight (with Battery): 2.3kg (5lb)

P-701 (701Series Pads)

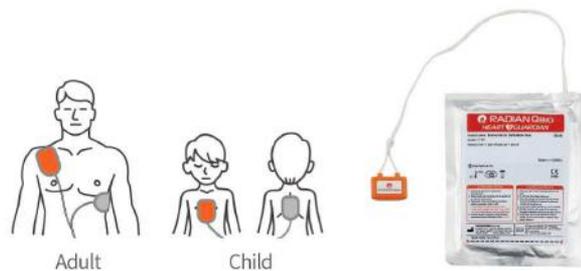

- Type: Pre-Gelled / Self-Adhesive / Disposable Non-Polarized Defibrillation Pads Non-Sterilization
- Shelf-Life: 2 years from Date of Manufacturing
- Size & Length: 110 x 120mm / 1.8m
- Pouch Size: 170 x 250mm

BT-701 (701Series Battery)



- Type: LiMnO₂(12V / DC6.4A) / Non-Rechargeable
- Weight: 285g
- Shelf-Life: 4years for Standby Mode (Can be reduced by temperature and Humidity condition during storage)

HR-B7 (701Series Semi-Rigid Carrying Case)



Closed

Open

Special Features

1 Energy



Adult : Min 150J~Max 250J(±15%)
Child : 50J(±15%)

2 Video Prompts



(4.3" High Resolution Color LCD)

3 Status Indicators

✓ Normal Operation



✓ Abnormal Operation



4 4 Languages support



5 Pre-connected Pads




Adult



Child

6 IP55 (Ingress Protection ratings)



7 Ambient Light Sensor



8 Communication Port

- ✓ USB Port 
- ✓ SD Card 
- ✓ Bluetooth 4.0 




- Semi-Automated External Defibrillator Biphasic Truncated Exponential (BTE) WaveForms
- Automatic Analysis ECG within 8 seconds
- Shows Patient Real Time ECG
- Adult / Child mode by Slide Switch
- Battery Power and Capacity 12V/6.4A (200Shocks or 10hrs continuous operation)
- Self-Tests(Daily, Weekly, Monthly) : Battery Status, Power Status, System Check Record
- Quick Power on & Charging the battery
- Review the data by PC Software (Option)

Accessories



Semi-Rigid Carrying Case

Pads

Battery

Manual